


# Stellisept<sup>®</sup> med gloves

## Antiseptic cleansing gloves



**Antiseptic body cleansing with comprehensive effect, e.g. against MRSA / ORSA and VRE**

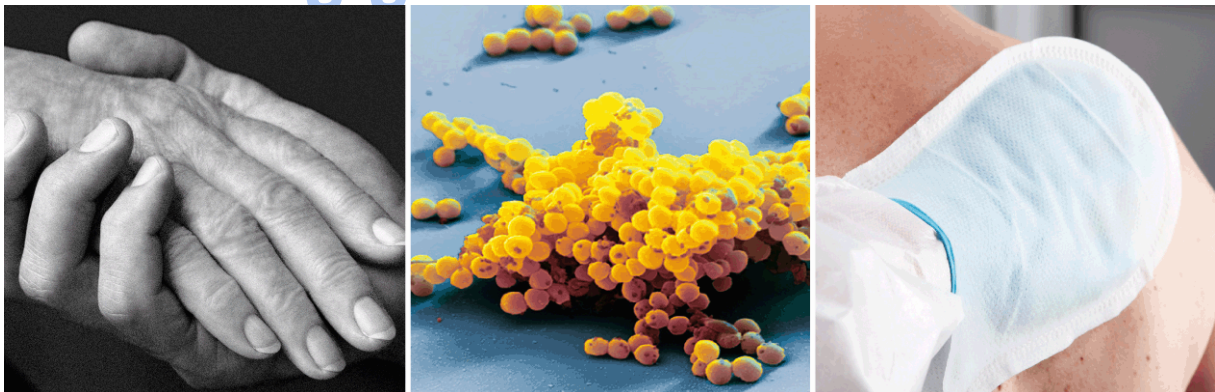

## Product properties

- comprehensive activity against national epidemic MRSA strains and clinical MRSA isolates
- cleansing gloves
- no rinse necessary after application
- skin nurturing properties
- colourant- and fragrance-free
- very good skin and mucous membrane tolerability
- can be heated in the microwave at 750 watt for 30 seconds

## Characteristics

Stellisept® med gloves are pleasantly soft PET gloves saturated with a mild, antiseptic cleansing solution. To make the whole-body wash as pleasant as possible, the gloves can be heated in the microwave (and can be reheated once) for approx. 30 seconds at 750 watt – or cooled for a refreshing treatment during summer.

## Composition

Aqua, Propyl Alcohol, Panthenol, PEG-40 Hydrogenated Castor Oil, Sodium Citrate, Didecylidimonium Chloride.

## Microbiological activity

- bactericidal (incl. MRSA / ORSA, VRE).

## Areas of application

Stellisept® med gloves are mild antiseptic cleansing gloves with replenishing properties and a broad spectrum of antiseptic effect. The gloves are suitable for the MRSA eradication in accordance with the RKI recommendations. The effect of Stellisept® med gloves was tested according to European standards. It meets the requirements of EN 1499 and EN 13727, and provides a comprehensive effect against national epidemic MRSA strains and clinical MRSA isolates. Through carefully selected ingredients, the cleansing gloves are very skin-friendly – a rinse with water after application is not necessary.

Stellisept® med gloves are suitable for the following areas of application:

### ■ antiseptic whole-body washing in case of MRSA / ORSA, VRE

Stellisept® med gloves possess both a proven comprehensive activity against MRSA and an excellent dermatological tolerability. The efficacy of Stellisept® med gloves against six national epidemic MRSA strains and numerous clinical isolates from different German regions was proven in tests analogous to the European standard EN 13727. Within 60 seconds, Stellisept® med gloves (saturation solution, 50%) were sufficiently active against all tested strains.

## Directions for use


One flow-pack contains 10 Stellisept® med gloves allowing you to use of one fresh glove per body region (Regions 1-10) to effectively prevent the spread of dirt or microorganisms.

Stellisept® med gloves are ready to use and must not be mixed with other soaps or solutions. Remove soap or wash lotion residues with water prior to the use of Stellisept® med gloves.

Long-term catheters can be also cleaned at the catheter entrance with a fresh Stellisept® med glove.

No rinse necessary after application.

- antiseptic whole-body washing of immobile patients in case of MRSA/ORSA

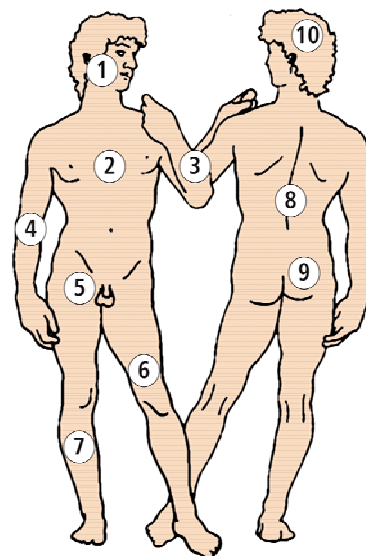

### Glove 1)

Clean the complete face with Stellisept® med gloves. The eyes need to be closed during this procedure. Treat the areas around the eyes with care, as the solution must not get into the eyes. When cleaning the ears, make sure the solution does not run into the ear.

### Glove 2)

For cleaning the torso, use Stellisept® med gloves and wash from the neck downwards.

### Gloves 3 - 4)

Wash arms from the fingertips to the shoulders. Make sure the axilla and gaps between fingers are sufficiently moistened with the preparation.

### Glove 5)

Use one glove for cleansing the genital area.

### Gloves 6 - 7)

Wash the legs from the hips downwards. Make sure the gaps between toes are sufficiently moistened with the preparation.

## Glove 8)

For cleaning the back, use Stellisept® med gloves and wash from the neck downwards.

## Glove 9)

Use one glove for cleansing the backside of the patient.

## Glove 10)

Stellisept® med gloves are also suitable for washing the hair. A subsequent rinse is not necessary.

No rinse necessary after application.

## Dosage

Depending on the size of the patient, use one flow-pack per patient and day for the MRSA eradication.

	Exposure time
Whole-body wash in case of MRSA / ORSA, VRE	60 sec.

## Safety instructions

Avoid contact with the eyes. In case of contact with eyes, rinse thoroughly with plenty of water, also under the eyelids. Seek medical advice in case the eyes remain irritated.

Only use undamaged packages. If warming up the gloves use them only lukewarm. Keep out of the reach of children. Protect against direct sunlight. Store product at temperatures below 30 °C.

## Chemical-physical data

Appearance	Tissue saturated with clear liquid
Smell	soapy, fresh
pH value (conc.)	5.5 – 6.5
Density (20 °C)	0.991 g/cm <sup>3</sup>

## Presentation

Flow-Pack contains 10 gloves

Note: The recommendations regarding our preparations are based on scientific tests and are given in good faith. More detailed recommendations, e.g. regarding material compatibility, are only possible in particular cases. Our recommendations are without obligation and do not constitute a warranty. They do not preclude a company's own testing for the intended purposes and processes. In this respect we cannot accept any liability. This complies with our general conditions of sale and supply.

08/11

